
# **O**Sleepnet<sup>®</sup>



featuring AIR gel with Advanced Cushion Technology

The Veraseal<sup>®</sup> 2 disposable mask brings AIR°gel technology to the acute care environment for a giant step forward in patient comfort. Three color coded elbow designs make choosing the right mask for CPAP/Bi-level, single or dual limb circuit NIV systems simple. The Veraseal 2 fits quickly and easily while enhancing patient comfort.

## Your Total Cost Solution

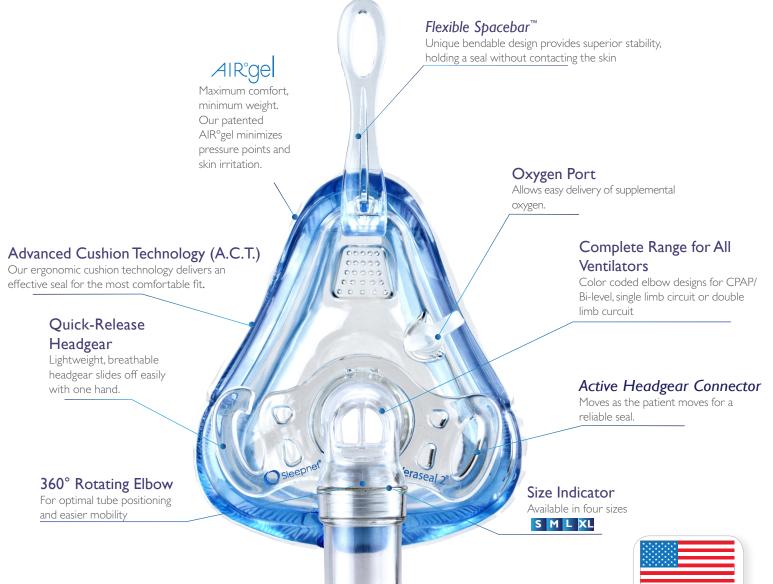

Made in the USA of US and imported parts

### What makes our Veraseal 2 Mask Advanced CushionTechnology(A.C.T.) different?

- Designed with minimal intentional leaks and elbows that work with CPAP/Bi-level and single or dual limb circuits.
- *Flexible Spacebar* is designed to provide superior stability without contacting the skin.
- Active Headgear Connector moves when the patient moves to help maintain an effective seal.
- Available in four sizes with simple and breathable quick-release headgear.



#### A full line of disposable masks to meet your NIV needs, color coded for easy identification



Vented Mask with an anti-asphyxia valve for use with CPAP or bi-level devices



Nonvented Mask with an anti-asphyxia valve for use with single limb circuits



Nonvented Mask for use with dual limb circuits

#### Veraseal 2 Masks with EZfit Disposable Headgear Specifications:

Product Number	Case Quantity	Style	Leak Port	Entrianment Valve	Size	Mask Weight	Mask Dimensions
50827	30	Vented	$\checkmark$	$\checkmark$	Small	84.3g	3.85''D × 3.90''W × 6.62''H
50828	30	Vented	$\checkmark$	$\checkmark$	Medium	89.2g	3.95''D × 4.21''W × 6.75''H
50829	30	Vented	$\checkmark$	$\checkmark$	Large	93.0g	4.05''D x 4.36''W x 6.75''H
50830	30	Vented		$\checkmark$	Extra Large	95.0 g	4.05''D x 4.50''W x 6.75''H
50933	30	Non-Vented	X	×	Small	84.3g	3.85''D x 3.90''W x 6.62''H
50934	30	Non-Vented	X	×	Medium	89.2g	3.95''D x 4.21''W x 6.75''H
50935	30	Non-Vented	X	×	Large	93.0g	4.05''D x 4.36''W x 6.75''H
50936	30	Non-Vented	X	×	Extra Large	95.0 g	4.05''D x 4.50''W x 6.75''H
50937	30	AAV Non-Vented	X	$\checkmark$	Small	84.3g	3.85''D x 3.90''W x 6.62''H
50938	30	AAV Non-Vented	X	$\checkmark$	Medium	89.2g	3.95''D x 4.21''W x 6.75''H
50939	30	AAV Non-Vented	Х	$\checkmark$	Large	93.0g	4.05''D x 4.36''W x 6.75''H
50940	30	AAV Non-Vented	X	$\checkmark$	Extra Large	95.0 g	4.05''D x 4.50''W x 6.75''H





# MiniMe<sup>®</sup>2

Sleepnet also offers pediatric masks for our smallest patients



# Contact us today.

To place your order, call an authorized distributor.

Call I.800.742.3646

Visit www.sleepnetmasks.com to learn more.

Sleepnet Corporation, 5 Merrill Industrial Drive, Hampton, New Hampshire 03842 USA

